

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USING

DANGER

To reduce the risk of electrocution:

- DO NOT use the unit while bathing.
- ALWAYS unplug this product immediately after using.
- DO NOT place or store the unit where it can fall or be pulled into a tub, sink, or toilet.
- DO NOT reach for a unit that has fallen into water; unplug immediately. Do not reach into the water under any circumstances. Have the unit checked by an Authorized Service Center after it has been in water.
- The unit is designed for use in a 120 volt AC, 60 Hz outlet. Never force the plug into the outlet.
- DO NOT submerge the unit in water or any other liquid.

- Never operate the unit if it has a damaged plug, if it is not working properly, if it has been dropped or damaged, or if it has been immersed in water. Refer to an Authorized Service Center for repair.
- DO NOT excessively bend or pull the cord. Unplug the cord by grasping the plug – NOT THE CORD.
- Keep the cord away from heated surfaces.
- Close supervision is necessary when this product is used by, on, or near children or disabled persons.
- Use this product only for its intended use as described in this manual. Use only attachments recommended by the manufacturer.
- Never use while drowsy.

WARNING

To reduce the risk of burns, electrocution, fire, or injury:

SAVE THESE INSTRUCTIONS

OPERATING INSTRUCTIONS

- This appliance has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit into a polarized outlet only one way. If the plug does not fit fully into the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature.

MEDICAL WARNING For your safety:

- If you have had oral or periodontal surgery in the previous twelve months, contact your dentist prior to use.
- When you start any new oral hygiene program, your gums may initially bleed slightly. Contact your dentist if excessive bleeding occurs or if bleeding continues after two weeks of use.
- Be sure that only the water jet nozzle touches your teeth and gums. Damage to teeth and gums may occur if any other part of the unit touches them.
- This product is designed for oral cleaning only. DO NOT use the unit for any other purpose.
- As with any other new oral hygiene program, individuals with heart problems should check with their physician prior to using the Interplak® Dental Water Jet, and should always exercise particular care with any oral hygiene device.
- The Interplak® Dental Water Jet is a personal care unit and is intended for household use only.

WARNING

Changes or modifications to this unit not expressly approved by Conair could void the user's authority to operate the equipment and void limited warranty. Servicing should always be referred to an Authorized Service Center.

Instructions For Use

Remove the water reservoir (2) and fill with lukewarm water. Reinstall the reservoir firmly onto the unit (1). Prior to initial use, the unit should be rinsed by letting a full reservoir of water pass through it. If the unit is to be used with a mouthwash solution, it should be rinsed with clean water after each use to remove any traces of mouthwash. If the unit has not been used for a long period of time, for reasons of hygiene it is advisable to follow the directions for initial use, letting a full reservoir of water pass through. To remove trapped air upon initial operation, turn the rotary flow control (7) all the way to the right and then to the left, three times.

Prior to filling the water reservoir, make sure water tip is firmly placed on and in the stand. The Water Jet is supplied with 5 attachments. Remove the tip support (3) and install selected tip:

Standard Tip (2pcs): removes trapped particles from between teeth.

Subgingival Tip (1pc): uses a soft rubber tip to massage the gums and safely deliver solution into the gingival pockets below the gumline.

Brush Head Tip (1pc): loosens debris and larger food particles before flushing them away – ideal for all dental work, including fixed and removable appliances.

Tongue Cleaner (1pc): place the tongue cleaner at farthest end of tongue and move it forward. Repeat to cover entire tongue. After switching on (6) the unit, the rotary flow control (7) can be used to regulate the water pressure. You should begin with the lowest setting so that your gums can become accustomed to the water jet.

Hold the fingertip on/off control (water-stop button) (4) down and, after switching on the unit, place the tip in your mouth and lean over the sink. When you release the water-stop button, the cleaning and massage process is begun. The direction of the water jet in your mouth can be controlled by twisting the tip to the desired position on the handle. Keep your mouth slightly open and hold your elbow at mouth level so that the water can run back into the sink. Direct the water jet toward your teeth, especially between your teeth, in order to rinse away stubborn remaining food.

This also massages your gums. Do not point the water jet directly at your gums; point between your teeth to prevent damage to delicate gum areas.

Limited 2-Year Warranty

Conair will repair or replace (at our option) your unit free of charge for 24 months from the date of purchase if the appliance is defective in workmanship or materials.

To obtain service under this warranty, return the defective product to the service centre listed below, together with your purchase receipt and \$5.50 for postage and handling. In the absence of a purchase receipt, the warranty period shall be 24 months from the date of manufacture.

ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE 24 MONTH DURATION OF THIS WRITTEN, LIMITED WARRANTY.

IN NO EVENT SHALL CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER.

This warranty gives you specific legal rights, and you may also have other rights, which vary from province to province.

© 2015 Conair Consumer Products ULC

SERVICE CENTRE:

Conair Consumer Products ULC
100 Conair Parkway, Woodbridge, Ontario, L4H 0L2.
1-855-926-6247

www.conaircanada.ca
Consumer Call Centre E-mail :
Consumer_Canada@Conair.com

15CN018049

IB-14087

INTERPLAK®

by CONAIR®

DENTAL WATER JET SYSTEM



Have a Question?

Please do not return this product to the retailer.

Call Us First!

Our customer service and product experts are ready to answer

ALL your questions.

Please call our TOLL-FREE customer service number at: 1-855-926-6247

or visit us online at www.conaircanada.ca


Instruction Manual • Model WJ7XBRC

For your safety and continued enjoyment of this product, always read the instruction book carefully before using.

www.interplak.com

The Interplak difference

The Interplak® Dental Water Jet was specifically designed to help you maintain a high level of oral hygiene and improved oral health quickly and easily at home. It flushes food debris from those hard-to-reach places around implants, orthodontic braces and crowns. It is gentle and effective on all orthodontic, crown, and bridge work.

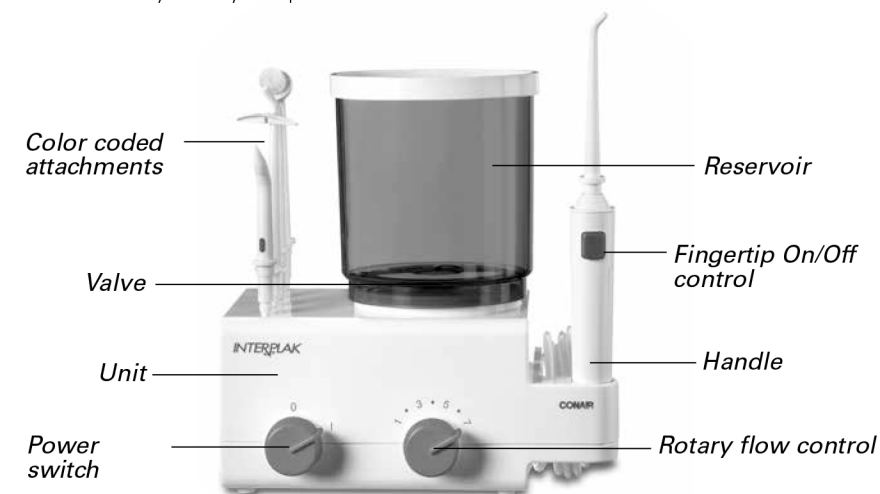

Cleaning and Care

If, in the course of use, the valve (8) on the water tank begins to leak, this can be taken care of by cleaning the valve seat. To do this, push out the valve body using the nozzle tip, clean the valve seat and replace the valve.



- 1 Unit
- 2 Reservoir
- 3 Handle
- 4 Fingertip on/off control
- 5 Color coded attachments

- 6 Power switch
- 7 Rotary flow control
- 8 Valve