



Congratulations on your purchase of the Interplak® Rotary Rechargeable Power Plaque Remover

Get ready to experience a brand new level of oral clean. State-of-the-art engineering has created a brush head with specially designed bristles that perform specific cleaning functions.

- Feathered subgingival bristles to clean under the gumline
- Contoured polishing tufts to smooth the tooth surface

With the push of one button, you get safe and effective three-way cleaning. Happy, healthy brushing, from Interplak!

Advanced Brush Head Design

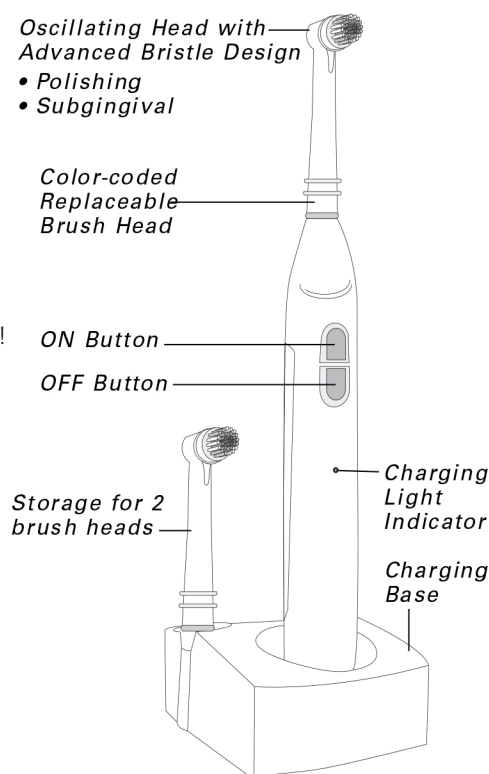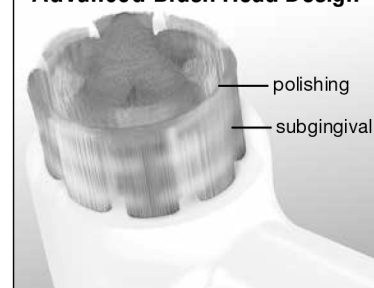

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

DANGER
To reduce the risk of electrocution:

- **READ ALL INSTRUCTIONS BEFORE USING**
- DO NOT handle the plug with wet hands.
- DO NOT place or store the charger base where it can fall or be pulled into a tub, sink, or toilet.
- DO NOT reach for a charger base that has fallen into water; unplug immediately. Do not reach into the water under any circumstances. Have the unit checked by an Authorized Service Center after it has been in water.
- The unit is designed for use in a 120 VAC, 60 Hz outlet. Never force the plug into the outlet.
- DO NOT submerge the charger base in water or any other liquid.

WARNING
To reduce the risk of burns, electrocution, fire, or injury:

- Close supervision is necessary when this product is used by, on, or near children or disabled persons.
- Never use while drowsy or while bathing.
- Keep the supply cord away from heated surfaces.
- Never operate the unit if it has a damaged plug, if it is not working properly, if it has been dropped or damaged, or if it has been immersed in water. Refer to an Authorized Service Center for repair.
- If the cord is damaged, discard the unit – the cord cannot be replaced.
- DO NOT excessively bend or pull the cord. Unplug the cord by grasping the plug – NOT THE CORD.
- Use this product only for its intended use as described in this manual. Only use attachments recommended by the manufacturer.
- This product is designed for teeth cleaning only – not for any other purpose.
- To prevent damaged bristles from breaking off while the unit is in use, replace your brush head if bristles are worn or bent.

SAVE THESE INSTRUCTIONS

MEDICAL WARNING
For your safety:

WARNING

Changes or modifications to this unit not expressly approved by Conair could void the user's limited warranty. Servicing should always be referred to an Authorized Service Center.

- Contact your dentist prior to use if you have had oral surgery in the previous twelve months.
- When you start any new oral hygiene program, your gums may initially bleed slightly. Contact your dentist if excessive bleeding occurs or if the bleeding continues after two weeks of use.
- Be sure that only the bristles touch your teeth and gums. Damage to teeth and gums may occur if any part of the unit touches them.
- This product is designed for oral cleaning only. DO NOT use the unit for any other purpose.
- As with any other new oral hygiene program, individuals with heart problems should check with their physician prior to using the RTG3XCSR, and should always exercise particular care with any oral hygiene device.
- Avoid using a brush with crushed or bent bristles. Damaged bristles may break off while the unit is in use. In order to prevent this, replace your brush head every three to four months, or sooner if the bristles become worn or bent.
- The RTG3XCSR is a personal care unit and is intended for household use only.

Initial Preparation


Before using, charge the power handle for 24 hours. To charge, press the off button on the power handle and place it in charger base. The indicator light will go on.

NOTE: Be sure the electrical outlet is functioning. You can test it with another appliance. Also, check the electrical outlet reset button and/or wall switch.

About twice a year, completely run down and then recharge your power handle to maximize the life of the rechargeable battery.


Instructions for Use: 5 Simple Steps

1. **Attach the brush head.**
After the unit has been fully charged, grasp the brush head by the neck and place the power handle shaft inside the hollow end of the brush head; (A) press down and (B) turn clockwise until the brush head snaps into place. To remove, grasp the brush head by the neck (not the bristles), turn counterclockwise, and pull away from the handle. **NOTE:** It is recommended that you replace your brush head every three months, or sooner if the bristles become worn or bent.
2. **Moisten the bristles.**
Hold the brush head under running water for a few seconds.
3. **Apply toothpaste sparingly and put brush head into mouth.**
4. **Turn the unit on.** Press the ON button.
5. **Position brush head properly.**
Place bristles perpendicular to flat tooth surface and guide the unit slowly along the outside, then the inside, and finally the biting surfaces of the teeth. You do not have to press hard.

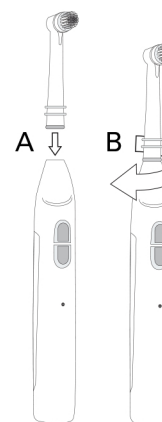

Cleaning and Maintenance of Unit


1. **After each use**
While toothbrush is on, hold the brush under running water for a few seconds.
2. **Turn the unit off**
Press the off button.
3. **Towel dry the power handle**
Place the handle onto the plugged-in charger base to maintain full power.
4. **Brush head maintenance**
Once a week or any time the power handle appears to be slowing down, turn power on and pour a small amount of mouthwash or water on the head.

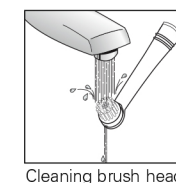

Battery Disposal

1. To remove the battery for disposal or recycling, insert a coin into one set of slots on the battery end cap and twist counterclockwise to release the battery cap.
2. Firmly tap the bottom of the handle on a hard surface in order to slide out battery.
3. Now the battery is exposed so you can separate the wire connectors. Carefully hold the wires on each side of the connector and gently cut the wires. The connector will separate and the batteries can be removed from the inner housing. Bring the battery to the nearest collection center for disposal or call the toll-free RBRC information line at 1-800-8-BATTERY. RBRC is a trademark of The Rechargeable Battery Recycling Corporation. DO NOT throw the batteries into the normal municipal waste system.

CONTAINS NICKEL METAL HYDRIDE BATTERY. MUST BE COLLECTED, RECYCLED OR DISPOSED OF PROPERLY. DO NOT DISPOSE OF IN THE MUNICIPAL WASTE.

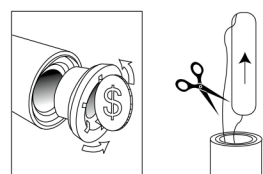

For more information on battery replacement for your power toothbrush, call 1-855-926-6247 (9am to 4:30pm Eastern Standard Time, Monday through Friday)

Every 6 months, completely run down and recharge your power handle to maximize the life of the rechargeable battery.



To help protect the environment, this product is powered by a rechargeable nickel metal hydride battery.

Under various state and local laws, it is illegal to dispose of this battery in our municipal waste stream. "The EPA certified RBRC Battery Recycling Seal on the nickel metal hydride (Ni-Mh) battery indicates Conair is voluntarily participating in an industry program to collect and recycle these batteries at the end of their useful life, when taken out of service in the United States or Canada. The RBRC program provides a convenient alternative to placing used Ni-Mh batteries in the trash or the municipal waste stream, which may be illegal in your area. Please call 1-800-9-BATTERY™ for information on Ni-Mh battery recycling and disposal bans/restrictions in your area. Conair's involvement in this program is part of our commitment to preserving our environment and conserving our natural resources."

*RBRC™ is a trademark of the Rechargeable Battery Recycling Corporation.

Power Brush Head Replacement

Dentists recommend replacing your toothbrush every 3-4 months and this applies to replacing the brush heads of your Interplak® Rechargeable Power Plaque Remover. You can find 2-Packs of Replacement Brush Heads, Model RBG3, at the same retail location where you purchased your power toothbrush. You can also find these replacement brush heads at most drug stores, discount and department stores, and food outlets in the oral care department. If you are having difficulty finding replacement brush heads, please call us at 1-855-926-6247 or visit us at www.interplak.com.

Product Features

Voltage: 120 VAC **Battery:** Nickel metal hydride rechargeable battery AAA type
Power: 3W (explosion may occur if replaced with other type battery)
Hertz: 60 Hz **Weight:** 230 grams

Limited 2-Year Warranty

Conair will repair or replace (at our option) your unit free of charge for 24 months from the date of purchase if the appliance is defective in workmanship or materials.

To obtain service under this warranty, return the defective product to the service centre listed below, together with your purchase receipt and \$5.50 for postage and handling. In the absence of a purchase receipt, the warranty period shall be 24 months from the date of manufacture.

ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE 24 MONTH DURATION OF THIS WRITTEN, LIMITED WARRANTY.

IN NO EVENT SHALL CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER.

This warranty gives you specific legal rights, and you may also have other rights, which vary from province to province.

© 2015 Conair Consumer Products ULC

SERVICE CENTRE:

Conair Consumer Products ULC
100 Conair Parkway, Woodbridge, Ontario, L4H 0L2.
1-855-926-6247

www.conaircanada.ca

Consumer Call Centre E-mail :
Consumer_Canada@Conair.com

15CN018442

IB-14174

INTERPLAK®

by CONAIR®

OPTICLEAN™ RECHARGEABLE
POWER PLAQUE REMOVER



Have a Question?

Please do not return this product to the retailer.

Call Us First!

Our customer service and product experts are ready to answer ALL your questions. Please call our TOLL-FREE customer service number at: 1-855-926-6247 or visit us online at www.conaircanada.ca

Instruction Manual • Model RTG3XCSRC

For your safety and continued enjoyment of this product, always read the instruction book carefully before using.

www.interplak.com