

MAINTENANCE

After about 12 paraffin treatments or when visible evidence of debris occurs, the unit should be cleaned and the paraffin wax should be replaced.

1. Unplug unit and remove lid and wax.

2. If wax is still hot, wait until it solidifies. Push on one end of the wax block to lift up.

DO NOT USE SHARP OBJECT. Remove wax and discard.

3. Clean inside of the unit with a light household cleaner and wipe dry.

4. Once unit is thoroughly clean and dry, place lid on unit and stow away. Do not add new paraffin wax until you are ready to use.

REPLACEMENT

CONAIR® Paraffin Spa Replacement Kit (Model #PWR5) is available through the Conair service center and retail stores carrying the PB8 unit. The kit includes 1.75 lbs of wax and 20 disposable liners.


WARNING:

IF YOU THINK THE POWER IS OFF WHEN THE SWITCH IS OFF, YOU'RE WRONG.

KEEP AWAY FROM WATER

Everybody knows that electricity and water are a dangerous combination. But did you know that an electrical appliance is still electrically alive even if the switch is off? If the plug is in, the power is on. So when you are not using your appliances, keep them unplugged.

DANGER — This appliance has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature.

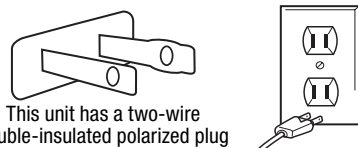

This unit has a two-wire double-insulated polarized plug

ALWAYS UNPLUG SMALL APPLIANCES

Public Service AD from Association of Home Appliance Manufacturers and Underwriters Laboratories, Inc.

LIMITED ONE YEAR WARRANTY

Conair will repair or replace (at our option) your unit free of charge for 12 months from the date of purchase if the appliance is defective in workmanship or materials.

To obtain service under this warranty, return the defective product to the service centre listed below, together with your purchase receipt. In the absence of a purchase receipt, the warranty period shall be 12 months from the date of manufacture.

ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE 12 MONTH DURATION OF THIS WRITTEN, LIMITED WARRANTY.

IN NO EVENT SHALL CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL, OR CON-SEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER.

This warranty gives you specific legal rights, and you may also have other rights, which vary from province to province.

SERVICE CENTRE

©2021 CONAIR CONSUMER PRODUCTS ULC

100 Conair Parkway

Woodbridge Ontario, L4H 0L2

For more information on any Conair

product call: 1-855-926-6247

Visit our Web site:

www.conaircanada.ca

Consumer Call Centre E-mail:

Consumer_Canada@Conair.com

21CN078106

renew · refresh · revive

true glow™

Thermal paraffin bath moisturizing system for healthy skin.

Soothing paraffin heat treatment



Have a Question?

Please do not return this product to the retailer.

Call Us First!

Our customer service and product experts are ready to answer **ALL** your questions.

Please call our TOLL-FREE customer service number at: 1-855-926-6247

or visit us online at www.conaircanada.ca

Instruction & Styling Guide

Model TGPB8C & TGPB8PRC

CONAIR®



IB-17231

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following: **READ ALL INSTRUCTIONS BEFORE USING THIS APPLIANCE.**

DANGER— To reduce the risk of electric shock:

- 1) Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
- 2) Do not reach for an appliance that has fallen into water. Unplug immediately.
- 3) Do not use while bathing or in a shower.
- 4) Do not place or store appliance where it can fall or be pulled into a tub or sink.
- 5) Do not place or drop into water or other liquid.

WARNING— To reduce the risk of burns, electrocution, fire, or injury to persons:

1. This product is for household use only.
2. This appliance should not be used by, on or near children or individuals with certain disabilities.
3. Unplug the appliance from outlet when not in use and before cleaning. Keep cord away from heated surfaces.
4. Do not operate this appliance with a damaged cord or plug or if the appliance malfunctions or has been dropped, dropped in water or damaged in any manner. Return appliance to nearest service facility for examination, electrical or mechanical adjustment or repair. No repairs should be attempted by the consumer.
5. Use appliance only for intended use as described in this instruction book. The use of attachments not recommended or sold by Conair may cause hazards.
6. Do not use outdoors or operate where aerosol (spray) products are being used, or where oxygen is being administered.
7. Never use while sleeping.
8. Do not overfill unit past MAX line inside unit. Use this appliance on level surface.
9. Do not heat paraffin wax in its plastic packaging and container.
10. Do not heat paraffin wax on stove or microwave. The paraffin unit is regulated at an even temperature to prevent scalding.
11. Do not operate unit without paraffin wax.
12. Do not pour melted paraffin wax from the unit.

13. Do not add water, lotions or any other additives to the wax.
14. Do not use if you have any skin irritations, rashes, cuts, burns, or warts.
15. Use only CONAIR® Paraffin Spa wax in this unit.
16. Do not soak or immerse unit in water. See MAINTENANCE for cleaning instructions.
17. Do not let the cord hang over the edge of the table or counter where it can be pulled on by children or tripped over.
18. Do not use extension cord with this appliance.
19. Never drop or insert any object into any opening.
20. This appliance has no serviceable parts.

SAVE THESE INSTRUCTIONS OPERATING INSTRUCTIONS

CAUTIONS

- This product should not be used in the following cases unless directed by a physician: benign and malignant tumors, hemorrhages, inflammations of the skin, phlebitis and thrombosis, open and fresh wounds, bruises, or broken skin.
- Consult your physician, should you have any questions regarding product's therapeutic use.
- In case of pregnancy, diabetes, or illness, consult your physician before using product.
- Do not use on any unexplained pain before consulting your physician.
- Use of product should be pleasant and comfortable. Should pain or discomfort result, discontinue use and consult your physician.
- Do not fall asleep while using this product.
- Keep away from children. This unit is not a toy. Children should not use it or play with it.

INSTRUCTIONS FOR USE

1. Place unit on a stable, level surface at a height that is convenient for treating the desired area.
2. Remove lid and wax container from inside the unit. Unwrap wax and discard all packaging.
3. Place the wax in the unit.
4. Place lid on unit. Lid should fit securely on unit.
5. Plug unit into electrical outlet (120V A.C./60 Hz). Switch on the unit. The indicator light will illuminate.

6. 1 lb of wax will take approximately 90 minutes to melt. Operate the appliance only with wax level below the MAX indicator located on the inside of the unit.

TREATMENT INSTRUCTIONS



1. Wash hands or feet thoroughly. Check to make sure that you have no cuts or visible patches of sensitive skin.
2. Remove any jewelry and roll up long sleeves.
3. If you are sensitive to temperature, test the wax before beginning. Simply dab a bit onto the inside of your wrist. It should be very warm to the touch, but not uncomfortable.
DO NOT DIP YOUR HAND IN WAX IF IT IS TOO HOT.
4. With fingers slightly apart, immerse hand slowly into wax in paraffin spa unit (**Fig.1**).
5. Immediately pull your hand out, being sure to maintain the open-finger position. Wait 3 to 5 seconds until the paraffin loses that “wet” look and then dip again. Try to immerse your hand to the same depth each time, so the lip of the wax “glove” thickens uniformly.
6. Repeat this process until 5 layers have been formed. Additional layers can be added to retain heat longer.
7. Insert your hand in one of the plastic liners provided.
8. Sit down and relax for 10-15 minutes, while the warm paraffin wax continues to release heat to soothe and soften the hand.
9. Remove plastic liner once you no longer feel any heat. Wiggle your fingers slightly to loosen the wax. Then peel off the paraffin starting from the wrist and pulling it downward like a glove, into the plastic liner (**Fig.2**).
10. Discard the plastic liner and wax. If only one person intends to use this unit, you may wish to recycle the wax by returning it to the unit. See replacement section on how to order replacement wax blocks and liners.
11. Return lid to unit, switch the unit to “0” and unplug. Allow unit and wax to cool before stowing away.